



ALPHA® LINERS

Bring comfort to your patients

- Unique interface materials
- Switch from gel to silicone within the same socket
- Customisable to challenging residual limb shapes



Alpha® Liners - Bring comfort to your patients

Use this chart to build your item number by going through the choices from top to bottom

AK

Select AK Symmetrical Liners

H440-1494 T440-1694-1 S440-1694-1

Type	Hybrid	SmartTemp	Silicone
	H	T	S
	Cushion	Locking Accordion Umbrella	
	440-	442-	

Gel profile	AK Symmetrical
Fabric	Hybrid Select Grey 4
Fabric	SmartTemp or Silicone Select Taupe 6

Gel thickness	AK Symmetrical 9mm distal 9					
Size Hybrid	4 M ₁	6 L	8 L ₊	7 XL		
Size SmartTemp Silicone	4-1 M ₁	4-2 M ₂	6-1 L ₁	6-2 L ₂	8-1 L ₊	7-1 XL
Dimensions	30-54 20-24	33-54 24-29	33-64 28-32	35-64 32-36	41-80 32-39	46-90 37-45

BK

Select Progressive Liners

H352-6394 T352-6594-1 S352-6594-1

Type	Hybrid	SmartTemp	Silicone
	H	T	S
	Cushion	Locking Accordion Umbrella	
	350-	352-	

Gel profile	Progressive
Fabric	Hybrid Select Grey / Light Grey 3
Fabric	SmartTemp or Silicone Select Taupe / Light Grey 5

Gel thickness	BK Progressive 2.5mm proximal 2.5mm posterior 6mm tibial crest 9mm distal 9									
Size Hybrid Classic	3 S	0 M	4 M ₁	6 L	8 L ₊	7 XL				
Size SmartTemp Silicone	3-1 S ₁	3-2 S ₂	0-1 M ₁	0-2 M ₂	4-1 M ₁	4-2 M ₂	6-1 L ₁	6-2 L ₂	8-1 L ₊	7-1 XL
Dimensions	19-32 14-22	21-32 16-23	26-44 20-24	29-44 24-29	30-54 20-24	33-54 24-29	33-64 28-32	35-64 32-36	41-80 32-39	46-90 37-45

Select Uniform Liners

H350-5364 T350-5360

Type	Hybrid	SmartTemp
	H	T
	Cushion	Locking
	350-	351-

Gel profile	Uniform
Fabric	Hybrid Select Grey / Light Grey 3
Fabric	Smart Temp Select Taupe / Light Grey 3

Gel thickness	BK Uniform 6mm 6					
Size Hybrid Classic	3 S	0 M	4 M ₁	6 L	8 L ₊	7 XL
Size SmartTemp Silicone	3 S	0 M	4 M ₁	6 L	8 L ₊	7 XL
Dimensions	21-32 15-22	23-36 18-26	33-50 20-28	36-55 25-33	40-60 28-40	43-65 33-45

Classic Original 6mm Uniform Liners

ALC-5064 ALL-5064

Type	Classic Cushion	Classic Locking
	ALC-	ALL-

Gel profile	Uniform
Fabric	Hybrid Select Grey / Light Grey 3
Fabric	Classic Original Buff 1

Gel thickness	BK Uniform 6mm 6					
Size Hybrid Classic	3 S	0 M	4 M ₁	6 L	8 L ₊	7 XL
Size SmartTemp Silicone	3 S	0 M	4 M ₁	6 L	8 L ₊	7 XL
Dimensions	21-32 15-22	23-36 18-26	33-50 20-28	36-55 25-33	40-60 28-40	43-65 33-45

Choose your material

WillowWood[®] ALPHA[®] CLASSIC

Alpha Classic gel consists of a skin-friendly blend of thermoplastic elastomer (TPE) with mineral oil and vitamin E, offering ultimate comfort. Its flow characteristics optimise pressure distribution, and its low rebound maintains its adapted shape over time, making this the perfect option for low to medium activity users.



FABRIC

Original

GEL PROFILE



Uniform

WillowWood[®] ALPHA[®] HYBRID

Designed to provide a more robust solution, the Hybrid gel is a thermoplastic elastomer (TPE), bridging Classic and Silicone thanks to its higher durometer. It is a perfect solution for users looking for a combination of comfort and durability.



Select



Uniform



Progressive



Symmetrical AK

WillowWood[®] ALPHA[®] SmartTemp

Combining silicone with the heat management technology of Outlast, the revolutionary SmartTemp gel regulates heat by responding to fluctuations in skin temperature. The silicone gel absorbs and stores heat from the limb, enhancing comfort and delaying the onset of sweating. As the body cools, the stored heat is released, maintaining a stable and comfortable socket environment.



Select



Uniform



Progressive



Symmetrical AK

WillowWood[®] ALPHA[®] SILICONE

Alpha Silicone utilises a platinum-cured medical grade silicone designed for comfort and performance, with an extremely high rebound for high levels of durability. A unique blend of vitamin E and mineral oils promotes skin repair and protection, providing a non-greasy and non-tacky interface for the user.



Select



Progressive




Symmetrical AK


GEL PROFILE




Uniform

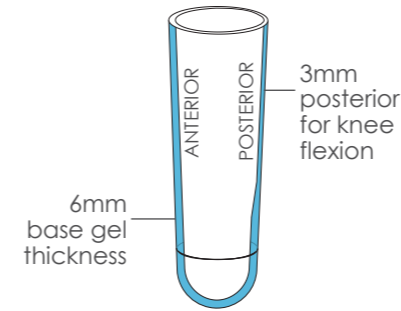

Progressive

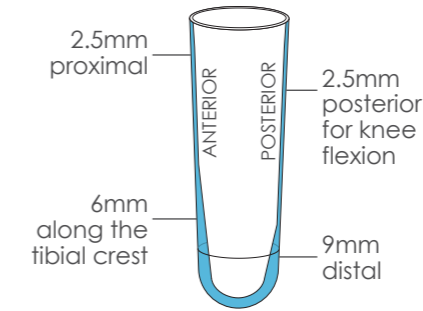

AK Symmetrical

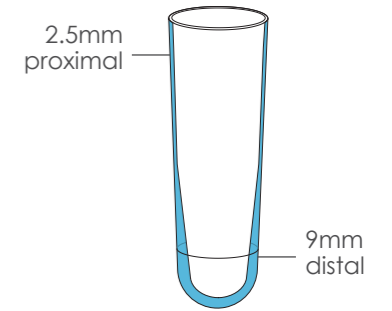

FABRICS

Select

Alpha Select Liners are constructed with a durable one-way stretch fabric. This helps to control pistoning without a distal matrix, whilst allowing radial stretch and maintaining elongation. BK liners feature a flexible Select panel, allowing easy bending and flexing of the knee.

Original

The popular all-purpose Original fabric ensures reliable comfort for the user, with a smooth texture to allow for easy socket donning.

RETROFITTING

Switch made easy

Alpha Liners with the same profile can retrofit into the same socket. This means patients don't need to switch socket to switch liners, making it easier to manage volume changes and accommodate a whole range of lifestyle activities.

ALPHA DESIGN LINERS

For patients with atypical residual limb contours, Alpha DESIGN Liners offer a tailored solution, with custom-designed gel profiles for each patient, providing a comfortable fit which is unmatched by off-the-shelf liners. DESIGN Liners are customised on the inside, providing a more typical exterior shape to fit a more standard socket shape.




UMBRELLAS

Uniform

Uniform profiles come with a standard fixed-diameter umbrella, relative to the size of the liner.

Accordion

Progressive and symmetrical profiles feature an accordion umbrella. The umbrella smoothly conforms to the distal end shape, ensuring an optimal, comfortable fit.

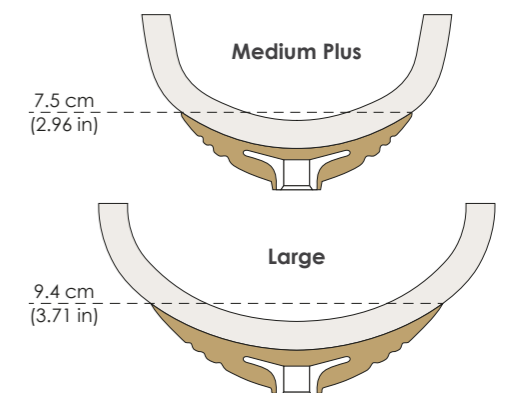

All Alpha Liners are available with or without a pin connection.